RēNEW
LIGHT THERAPY

7LS[™]
BEAUTY WELLNESS

PLEASE TAKE A
MOMENT NOW

REGISTER YOUR PRODUCT AT:

www.7LSbeauty.com

YOUR VALUABLE INPUT REGARDING THIS
PRODUCT WILL HELP US CREATE
THE PRODUCTS YOU WILL WANT
IN THE FUTURE.

INSTRUCTION MANUAL AND WARRANTY INFORMATION

1-YEAR LIMITED WARRANTY

FA7-1450J

READ ALL INSTRUCTIONS BEFORE USE.

SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE.

WARNINGS, PRECAUTIONS, CONTRAINDICATIONS:

FOLLOW THE INSTRUCTIONS OUTLINED IN THIS GUIDE CLOSELY TO ENSURE YOUR SAFETY AND FLAWLESS OPERATION:

- This product is only intended for use on the face.
- Discontinue using the product immediately if any abnormal situation, allergic reaction, or discomfort occurs during use.
- If you are under medical treatment, please consult with your doctor before use.
- Keep the device clean after using.
- DO NOT expose to eyes directly.
- Before operation, wear eye protection goggles.
- DO NOT use any detergent or liquid to clean the surface of this product.
- This product has a special lacquer. When cleaning, do not use any corrosive solution which may cause damage to the lacquer.
- Repair, removal, or modification of this product is strictly prohibited. Improper handling may result in short circuits, product damage, or even fire. Maintenance and warranty services will be void if unauthorized changes have been made to the product. Please contact Consumer Relations if the product has a problem.
- Be cautious of any edges, uneven surfaces, metal parts, accessories, or packaging to prevent possible injury or damage. Take note that this product should not be placed in your mouth.
- DO NOT use this product in high temperatures or high humidity levels (e.g., bathrooms). Ideally, the product should be used in a location free of excessive dust and not exposed to direct sunlight.
- DO NOT insert foreign objects or liquids into the product as this could damage the product.
- Refrain from dropping the product or bumping it against other objects. Avoid damage, collision, and excessive pressure by always handling the product with care.
- DO NOT place this product near strong magnetic or electric fields which may cause product failure.
- DO NOT use this product in hazardous locations.
- **CAUTION: The product may explode if disposed in a fire.**
- Due to hygienic reasons, do not share this product with other people.

CHILDREN

Keep out of reach of children. Do not allow children to use or play with product or its accessories. Product and accessories may contain small parts that could become detached and create a choking hazard.

DISPOSAL OF ELECTRICAL AND ELECTRONIC EQUIPMENT


This symbol indicates that all electrical and electronic equipment included with this product should not be treated as household waste. Instead, it should be left at the appropriate collection point for the recycling of electronic equipment. By ensuring this product is disposed of correctly, you are helping to prevent negative consequences for the environment and human health. To recycle this product, please contact your local city office, environmental protection agency, household waste disposal service or where you purchased the product to learn more about appropriate recycling methods.

BATTERY

For best battery performance, please allow battery to deplete 100% the first 3 uses.

A rechargeable battery can be charged and discharged over 100 times. However, it will eventually wear out. This is not a defect and corresponds to normal everyday use. When the operating or standby time is noticeably shorter, it may be a sign that the battery has reached the end of its life span. Please do not make repairs, modifications, or installations.

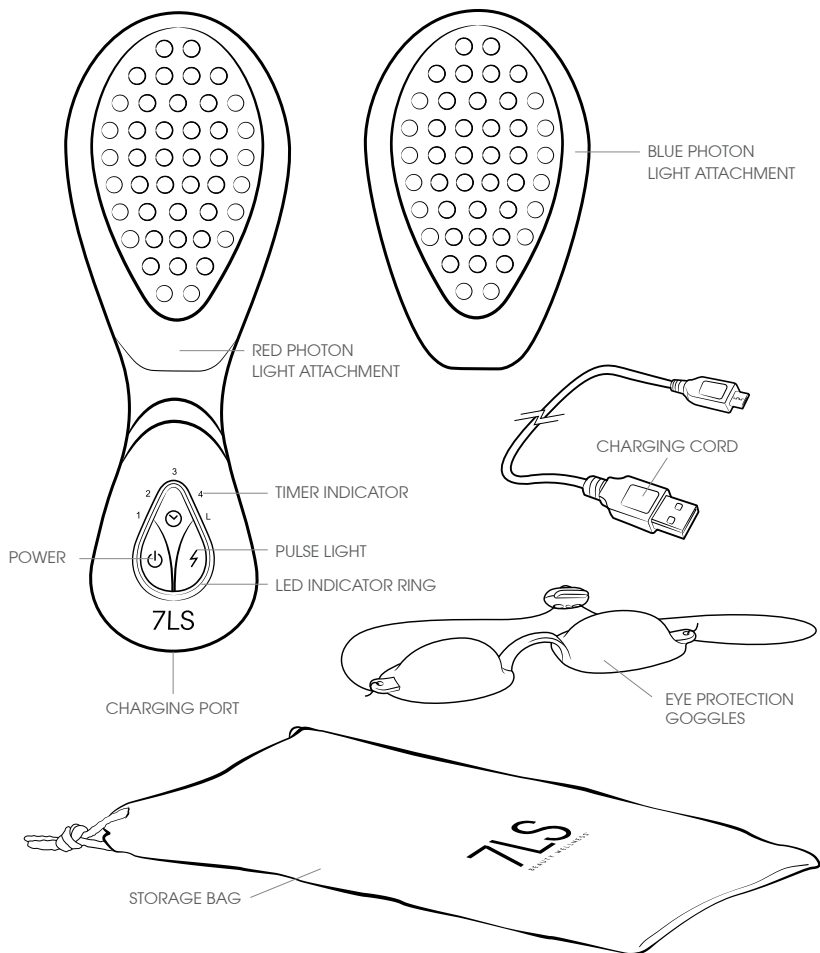UNIQUE FEATURES & SPECIFICATIONS

Breathe new life into your skincare routine with the introduction of light therapy. **RēNEW** boasts a noninvasive LED treatment featuring complexion-clarifying elements to minimize wrinkles, smooth skin tone and texture, and reduce blemishes for in-office-worthy results—with no down time.

INTERCHANGEABLE ATTACHMENTS

For different treatment types

RED LIGHT TREATMENT

Helps activate skin cells by stimulating collagen

BLUE LIGHT TREATMENT

Helps balance oil production and reduce inflammation

TIMED TREATMENT MODES

Provides timed treatment for different skin care needs

RECHARGEABLE

Charging cord

HOW TO USE

PRIOR TO USE

- Make sure the unit is fully charged.
- Make sure the treatment head of choice is attached.

CHARGING

- Connect charging cord to charging port.
- Connect USB to a power supply.

BATTERY INDICATOR


Low battery: Blinking red light on LED indicator ring

During charging: Red light on LED indicator ring

Battery full: No light on LED indicator ring

HOW TO CHANGE THE TREATMENT HEAD

- Hold the unit handle in one hand while pinching the top of the connected treatment attachment in the other. Slide the attachment away from the unit (you will feel a slight magnetic resistance).
- Continue holding the unit handle in one hand and slide the newly chosen treatment attachment onto the top of the unit.

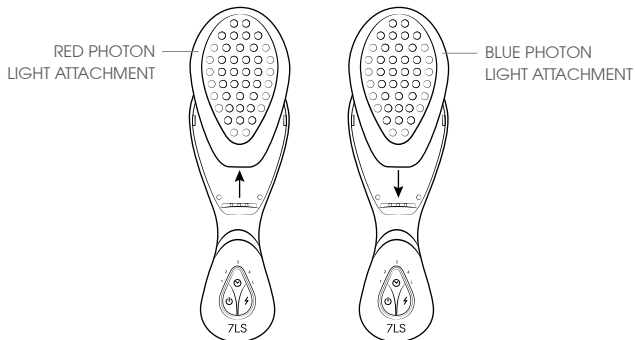HOW TO USE

RēNEW RED LIGHT TREATMENT


- Thoroughly wash and dry face.
- Make sure the red photon light is attached.
- Put on the eye protection. **CAUTION: Do not use device without wearing eye protection goggles. Damage to your eyes may occur.**
- Press **⏻** (blue LED is always on) to start. A beep will sound. Lights default to the last selected working mode.
- Hold the lights approximately 2 inches from the face and move slowly from chin to forehead.
- For best results, use this treatment for 3 minutes per treatment area, 2 times per week for 6 weeks.

RēNEW BLUE LIGHT TREATMENT

- Thoroughly wash and dry face.
- Make sure the blue photon light is attached.
- Put on the eye protection. **CAUTION: Do not use device without wearing eye protection goggles. Damage to your eyes may occur.**
- Press **⏻** (blue LED is always on) to start. A beep will sound. Lights default to the last selected working mode.
- Hold the lights approximately 2 inches from the face and move slowly from chin to forehead.
- For best results, use this treatment for 3 minutes per treatment area, 2 times per week for 6 weeks.

TREATMENT TIME SETTINGS

- In treatment mode, press ⌚, the first level blue LED will turn on. After 1 minute, the device will beep 3 times consecutively and turn off.
- In the first level position, press ⌚, the second level blue LED will turn on. After 3 minutes, the device will beep 3 times consecutively and turn off.
- In the second level position, press ⌚, the third level blue LED will turn on. After 5 minutes, the device will beep 3 times consecutively and turn off.
- In the third level position, press ⌚, the fourth level blue LED will turn on. After 12 minutes, the device will beep 3 times consecutively and turn off.
- In the fourth position, press ⌚, the L level blue LED will turn on. After 120 minutes, the device will beep for 3 times consecutively and turn off.
- In the L position, press ⌚ to enter standby mode. After 1 minute, the device will beep 3 times consecutively and turn off.

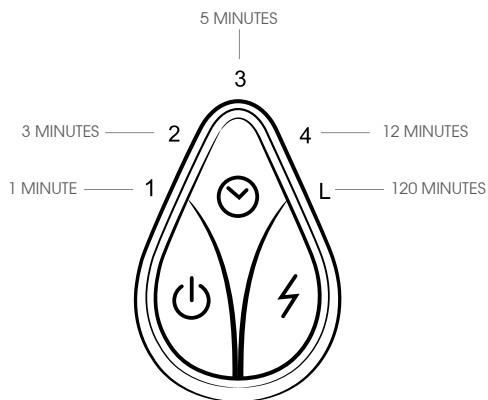RENEW PULSE LIGHT MODE

- In any operating mode, press  once (blue LED is on) to enter pulsed mode. Lights will blink.
- Press  (blue LED goes off) again for a beep sound, return to the working mode, light is on.

POWERING OFF

- In any operating mode, press , it will beep 3 times consecutively and then turn off.
- In any operating mode, after setting the timing, the device will beep 3 times consecutively after set time has elapsed.

CLEANING AND CARE

CAUTION: Before cleaning the unit, turn power off.

CAUTION: Use of detergents, liquids, or other cleaning products other than a soft, dry cloth may cause damage to the product.

Clean the LED surface with a dry cloth after use.

TROUBLESHOOTING

PROBLEM	POSSIBLE CAUSE	SOLUTION
Unit will not operate	• Unit is not charged • Power button not pressed • No power at outlet	• Charge the unit • Press and hold the power  button for 1 second to turn the power on • Check circuits, fuses. Try a different outlet

FAQ

WILL IT HURT?

Treatment should not be painful. Depending on the sensitivity of your skin, you may feel slight discomfort from the warmth of the LEDs. There could be some short-term redness or sensitivity after treatment. If treatment is uncomfortable, you may need to use a lower treatment setting.

CAN I USE THE PRODUCT IN COMBINATION WITH ANTIAGING CREAMS/SERUMS?

If your antiaging cream or serum contains AHA (or acid), then it should not be used on the same day as the treatment, but can be used on nontreatment days. Non-AHA (acid free) antiaging products can be used well before treatment, to allow absorption.

1-YEAR LIMITED WARRANTY

FKA Brands sells its products with the intent that they are free of defects in manufacture and workmanship for a period of 1 year from the date of original purchase, except as noted below. FKA Brands warrants that its products will be free of defects in material and workmanship under normal use and service. This warranty extends only to consumers and does not extend to Retailers.

To obtain warranty service on your FKA Brands product, contact a Consumer Relations representative for assistance. Please make sure to have the model number of the product available.

FKA Brands does not authorize anyone, including, but not limited to, Retailers, the subsequent consumer purchaser of the product from a Retailer or remote purchasers, to obligate FKA Brands in any way beyond the terms set forth herein. This warranty does not cover damage caused by misuse or abuse; accident; the attachment of any unauthorized accessory; alteration to the product; improper installation; unauthorized repairs or modifications; improper use of electrical/power supply; loss of power; dropped product; malfunction or damage of an operating part from failure to provide manufacturer's recommended maintenance; transportation damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of FKA Brands.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adoption to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. THERE SHALL BE NO OTHER WARRANTIES EXPRESS OR IMPLIED INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS OR ANY OTHER OBLIGATION ON THE PART OF THE COMPANY WITH RESPECT TO PRODUCTS COVERED BY THIS WARRANTY. FKA BRANDS SHALL HAVE NO LIABILITY FOR ANY INCIDENTAL, CONSEQUENTIAL, OR SPECIAL DAMAGES. IN NO EVENT SHALL THIS WARRANTY REQUIRE MORE THAN THE REPAIR OR REPLACEMENT OF ANY PART OR PARTS WHICH ARE FOUND TO BE DEFECTIVE WITHIN THE EFFECTIVE PERIOD OF THE WARRANTY. NO REFUNDS WILL BE GIVEN. IF REPLACEMENT PARTS FOR DEFECTIVE MATERIALS ARE NOT AVAILABLE, FKA BRANDS RESERVES THE RIGHT TO MAKE PRODUCT SUBSTITUTIONS IN LIEU OF REPAIR OR REPLACEMENT.

This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on internet auction sites and/or sales of such products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate as to any products or parts thereof which are repaired, replaced, altered, or modified, without the prior express and written consent of FKA Brands.

This warranty provides you with specific legal rights. You may have additional rights which may vary from country to country. Because of individual country regulations, some of the above limitations and exclusions may not apply to you. For more information regarding our product line in the USA, please visit: www.fkabrand.com. For Canada, please visit: www.fkabrand.ca.

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